

Bugbox

ANAEROBIC & MICROAEROPHILIC WORKSTATIONS GLOBAL (220:110V)

Bugbox anaerobic workstations are designed specifically to help microbiologists cope with rising workloads and provide the best primary isolation rates.

A NEW LEVEL OF VERSATILITY, EFFICIENCY AND POWER IN A COMPACT SPACE

Bugbox is easy to use. Its compact size meets the needs of even the smallest laboratory spaces.
Adjustable temperature and humidity provides a precisely controlled anaerobic environment that is optimal for cell growth, with no dry spots.



Experience Baker Ruskinn Anaerobic Workstations



Compared with approximately 20 anaerobic jars per week, Bugbox is economical with a lower cost per plate, more reliable, providing a stable atmosphere, and minimal maintenance.

DESIGNED TO PROTECT YOUR RESULTS

- The acrylic airtight chamber is flooded with anaerobic gas mix (H₂ in N₂) and O₂ is displaced.
- If any O₂ remains or is allowed to enter, it is "scavenged" by a palladium catalyst situated under the floor tray the O₂ reacts with the H₂ to form water.
- Interlock uses an N₂ purge, so when a user brings in plates through the interlock, no O₂ enters the main chamber
 inner and outer interlock doors cannot be opened simultaneously.
- Gloveless Ezee Sleeves[™] are purged using N₂ gas via foot pedals, so no O₂ enters the main chamber when the glove ports are opened.

ECONOMIC AND RELIABLE FOR LONG TERM SAVINGS

- Standard dual gas operation, low gas consumption and running costs.
- Lower cost per plate compared to anaerobic jars.
- Minimal maintenance and downtime.

ULTIMATE CONTROL FOR OPTIMUM CELL ENVIRONMENT


- Accurate temperature control from ambient + 5°C to 45°C.
- Accurate and automated humidity control, no dry spots.
- Palladium catalyst maintains anaerobic environment, plus anaerobic colorindicator strips verify anoxic conditions.
- Ezee Sleeve™ Direct Hand entry system allows access without disrupting the atmosphere within the chamber.

CONVENIENT & COMFORTABLE USER EXPERIENCE

- 1 QUICK AND EASY DIRECT ACCESS
- 2 SHORTEST INTERLOCK CYCLE TIME IN THE INDUSTRY
- 3 SINGLE PLATE ENTRY SYSTEM (SPES)
- **4** ENERGY-SAVING LIGHTING
- 5 HIGH-INTENSITY INSPECTION SPOT LAMP

MICROAEROPHILIC OPTIONS AVAILABLE



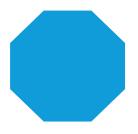

OUR ANAEROBIC
WORKSTATIONS ARE
DESIGNED SPECIFICALLY TO
HELP MICROBIOLOGISTS COPE
WITH RISING WORKLOADS
AND PROVIDE THE BEST
PRIMARY ISOLATION RATES.



BUGBOX SPECIFICATION SUMMARY

MODEL		BUGBOX/ BUGBOX M	BUGBOX PLUS
External Dimensions	Width	800 mm	830 mm
	Depth	660 mm	660 mm
	Height	650 mm	650 mm
Internal Dimensions	Width	500 mm	500 mm
	Depth	460 mm	460 mm
	Height	420 mm	420 mm
Maximum Capacity	90 mm Plates	270	234
Working Capacity	90 mm Plates	200	180
Interlock Dimensions	Width	100mm	150 mm
	Depth	100mm	230 mm
	Height	200mm	190 mm
Interlock Capacity	90 mm Plates	10	18
Interlock Time Cycle		15 sec	35 sec
Interlock Door Operation		Manual	Manual
Weight*		99 lbs / 143 lbs	121 lbs
Petri Dish Holders (Standard)		3	3

^{*}Bugbox M (143lbs) includes ICONIC™, the gas mixing system from Baker Ruskinn that gives Bugbox M unparalled oxygen control while using up to 40% less nitrogen.

STANDARD FEATURES


- Detox advanced carbon filtration system
- Ezee Sleeve $^{\mathrm{TM}}$ direct hand entry system
- Energy saving fluroescent illumination
- Inspection spot lamp
- Low gas alarm
- Automatic humidity control
- Palladium catalyst
- Anaerobic indicator strips
- Petri dish holders

OPTIONS & ACCESSORIES

- Vacuum line port
- Gas sample port
- Cable gland port
- Internal electrical outlet
- Gas tank regulators and filter modules
- Workstation stand
- External docking facility for anaerobic jars
- Power failure back-up system
- Data logging connection
- Single Plate Entry System (SPES)
- For faculative and microaerophiles
 Bugbox-M allows user defined control of O₂ and CO₂ (featuring ICONIC™ from Baker Ruskin)

ICONIC™ allows:

- O₂ from 0.0% to 23.0% in 0.1% increments
- CO₂ control from 0.1% to 30.0% in 0.1% increments
- O₂ sensor calibration with one touch
- Microaerophilic cycling, facilitating up to 4 different O₂ and CO₂ concentrations through a user-defined sequence of time

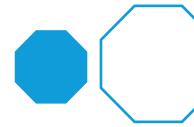

Get in touch today to hear how we can work with you.

For U.K. / global inquiries (except U.S./Canada): sales@ruskinn.com +44 (0) 1656 645988

Connect and Collaborate with Baker









Anaerobic & Microaerophilic Workstations

If you are looking for more robust processing power and capacity, the Baker Ruskinn Concept range of workstations are the perfect addition to any lab.







VERSION NO.001

www.bakerco.com

Designed in the U.K.